

Methods for Collecting, Detecting and Identifying Biological Agents in Environmental and Animal Samples (Briefing Charts)



Johnathan L. Kiel,
Eric A. Holwitt, Veronica K. Sorola,
Maomian Fan,
Yvette B. Gonzalez, Ishmael I. Rosas
and David F. Vela

711th Human Performance Wing
Human Effectiveness Directorate
US Air Force Research Laboratory

Report Documentation Page				Form Approved OMB No. 0704-0188	
Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.					
1. REPORT DATE 01 AUG 2009		2. REPORT TYPE		3. DATES COVERED	
4. TITLE AND SUBTITLE Methods for Collecting, Detecting & Identifying Biological Agents in Environmental and Animal Samples				5a. CONTRACT NUMBER In-House	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER 62202F	
6. AUTHOR(S) Johnathan Kiel; Eric Holwitt; Veronica Sorola; Maomian Fan; Yvette Gonzalez				5d. PROJECT NUMBER 7757	
				5e. TASK NUMBER P4	
				5f. WORK UNIT NUMBER 01	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) AFMC,AFRL,711 HPW, Human Effectiveness Directorate,Biosciences and Performance,Division, Biobehavior, Bioassessment, & Biosurveillance Branch,Brooks City-Base,TX,78235				8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release; distribution unlimited.					
13. SUPPLEMENTARY NOTES					
14. ABSTRACT N/A					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT	18. NUMBER OF PAGES 31	19a. NAME OF RESPONSIBLE PERSON
a. REPORT unclassified	b. ABSTRACT unclassified	c. THIS PAGE unclassified			

Outline

- Advantages of aptamers
- ALISA, Where we came from
- One step Quantum Dot De-quenching Assay
 - Why we need a double-stranded DNA aptamer
- Comparing **SELEX** to Aptamer Selection Express (**ASExpP**)
- Reagentless electronic sensors (RFIDs)
- Emerging disease agents and finding an unknown
 - Why we need a rapid technique for aptamer selection
- Summary

Advantages of Aptamers

- Aptamers are smaller than antibodies – ranging from 30 to 50 nucleotides
- Do not require either animals or tissue culture for production
- Can be synthesized chemically or by PCR
- Due to the nature of DNA, they are stable in harsh environments and do not require special storage conditions
- Offer additional chemistries and modalities for further stabilization (nuclease resistance) and assays

Where we came from: Tularemia in Houston: PCR and Immunoassays are not the last word

Berger, “Suspicious bacteria detected: Security monitors spot germ; terrorism discounted,” *The Houston (TX) Chronicle* 10 October 2003:A27

Francisella tularensis also discovered on Washington (DC) National Mall 24-25 Sept 2005; not reported until 1 Oct 2005

Laboratory Investigation (2006) 1-9
© 2006 USCAP, Inc. All rights reserved 0023-6837/06 \$30.00
www.laboratoryinvestigation.org



Technical Report

Anti-*Francisella tularensis* DNA aptamers detect tularemia antigen from different subspecies by Aptamer-Linked Immobilized Sorbent Assay

Jeevalatha Vivekananda and Johnathan L. Kiel

Air Force Research Laboratory, HEPC, Brooks City-Base, TX, USA

Aptamers are powerful candidates for molecular detection of targets due to their unique recognition properties. These affinity probes can be used to recognize and bind to their targets in the various types of assays that are currently used to detect and capture molecules of interest. They are short single-stranded (ss) oligonucleotides composed of DNA or RNA sequences that are selected *in vitro* based on their affinity and specificity for the target. Using combinatorial oligonucleotide libraries, we have selected ssDNA aptamers that bind to *Francisella tularensis* subspecies (*subsp.*) *japonica* bacterial antigen. *F. tularensis* is an intracellular, nonmotile, nonsporulating, Gram-negative bacterial pathogen that causes tularemia in man and animals. Just as antibodies have been used to detect specific targets in varying formats, it is possible that nucleic acid-binding species or aptamers could be used to specifically detect biomolecules. Aptamers offer advantages over antibody-based affinity molecules in production, regeneration and stability due to their unique chemical properties. We have successfully isolated a set of 25 unique DNA sequences that specifically bind to *F. tularensis* subspecies *japonica*. When tested in a sandwich Aptamer-Linked Immobilized Sorbent Assay (ALISA) and dot blot analysis, the aptamer cocktail exhibited specificity in its ability to bind only to tularemia bacterial antigen from subspecies *japonica*, *holarctica* (also known as *palaearctica*) and *tularensis* but not to *Bartonella henselae*. Moreover, there is no binding observed either to pure chicken albumin or chicken lysozyme. Thus, it appears that this novel antitularmia aptamer cocktail may find application as a detection reagent for a potential biological warfare agent like *F. tularensis*.

Laboratory Investigation advance online publication, 20 March 2006; doi:10.1038/labinvest.3700417

Keywords: DNA aptamers; *Francisella tularensis*; SELEX; ALISA; ELISA; dot blot; bioterrorism

Aptamers are single-stranded oligonucleotides with a length of tens of nucleotides, obtained by systemic evolution of ligands by exponential enrichment (SELEX) technology, exhibiting high affinity and specificity towards any given target molecule.^{1,2} These single-stranded (ss) nucleic acid molecules have highly defined tertiary structures, which allow them to form stable and specific complexes with a range of different targets including small molecules such as amino acids to highly complex proteins and whole viruses.³⁻⁷ For example, DNA-binding species have been selected that can interact with thrombin⁸ and RNA aptamers have been selected that recog-

nize a variety of cytokines.⁹ These specialized molecules are analogous to antibodies in specificity and affinity with an apparent advantage of being reproduced by chemical synthesis and more easily labeled with fluorescent or other reporters during their synthesis. Comparisons of various ligand-binding aptamers with proteins that bind to their targets have shown that both nucleic acids and proteins use similar strategies for the formation of well-defined binding patterns.^{10,11} Structural studies with aptamer-target complexes have demonstrated insights into molecular diversity associated with nucleic acid architecture and molecular recognition.¹² Aptamers frequently form complexes that have dissociation constants in the nanomolar range and can clearly distinguish between even closely related protein targets.^{13,14} The degree of molecular distinction achieved by aptamers may surpass that of antibodies with a remarkable diversity in structure and function.¹⁵

Correspondence: Dr J. Vivekananda, PhD, HEPC, Air Force Research Laboratory, 2486 Gillingham Dr, Bldg 175 E, Brooks City Base, TX 78235, USA.
E-mail: jeevalatha.vivekananda@brooksa.af.mil
Received 9 September 2005; revised and accepted 10 February 2006; published online 20 March 2006

Approved for public release; distribution unlimited

Current Methods for Tularemia Diagnosis

Culture in cysteine enriched medium with glycerol (dangerous for lab personnel)

1-10 Organisms can cause infection

Type A (most pathogenic)--glycerol catabolism positive (exceptions: *Ft mediasiatica* and *novicida*)

Type B (self-limiting)--glycerol catabolism negative

CDC PCR method

Smaller product--**Type A**

Larger product (insertion)--Type B

Immunoassays (ELISA and Agglutination)



For antibodies to *F. tularensis* in serum

Can be adapted to find bacterial antigen


Not type specific

The Houston cat after preliminary culture



Type B

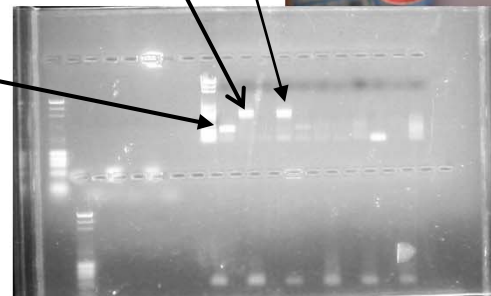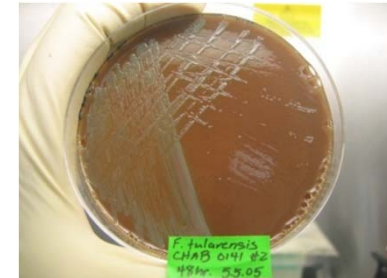

Anti - Francisella tularensis Polyclonal Antibody - ELISA

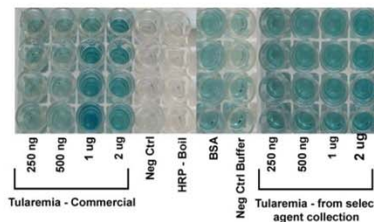

Anti-Francisella tularensis Aptamer ELISA

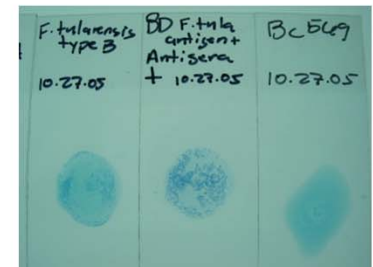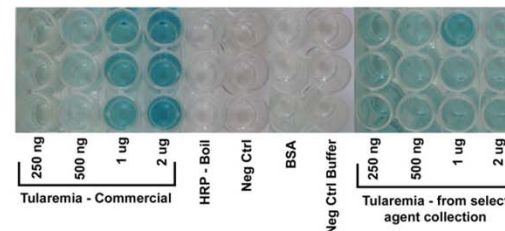

Approved for public release; distribution unlimited

Sensitivity of Aptamers for Detecting *Bacillus thuringiensis* Spores and *Francisella tularensis*

J Fluoresc (2007) 17:193–199
DOI 10.1007/s10895-007-0158-4

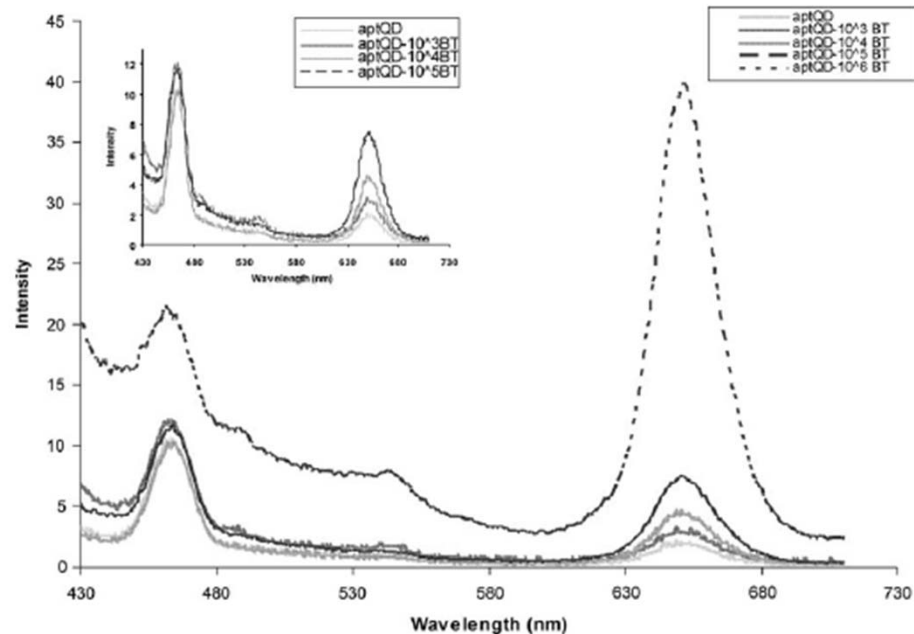


Fig. 3 Fluorescence Spectra of aptamer-QD bound to BT spores. Dashed line is the spectrum of aptamer-QD reacted with 10^6 CFU of BT (apt QD- 10^6 -BT). Other dilutions are shown in the legend area. An inset in the upper left hand corner represents the same data with the exclusion of the data for 10^6 CFU of aptamer-QD bound to BT spores

Laboratory Investigation advance online publication, 20 March 2006; doi:10.1038/labinvest.3700417

DNA aptamers for *F. tularensis* antigen
J Vivekananda and J.L. Kiel

9

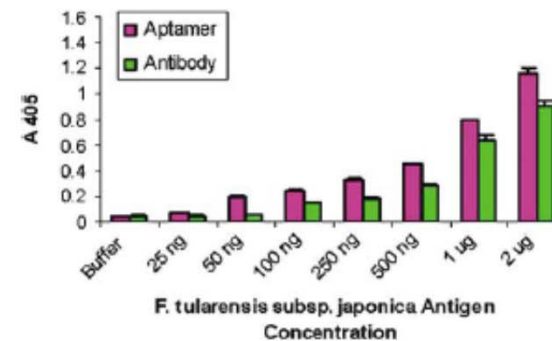


Figure 1 Sensitivity of anti-tularemia aptamer cocktail for *F. tularensis* subspecies *japonica* antigen and anti-tularemia antiserum as assessed by ALISA and ELISA. The assays were performed as described in 'Materials and methods'. The data are presented as OD at 405 nm vs antigen quantity. Averages of four replication measurements are shown in the figure.


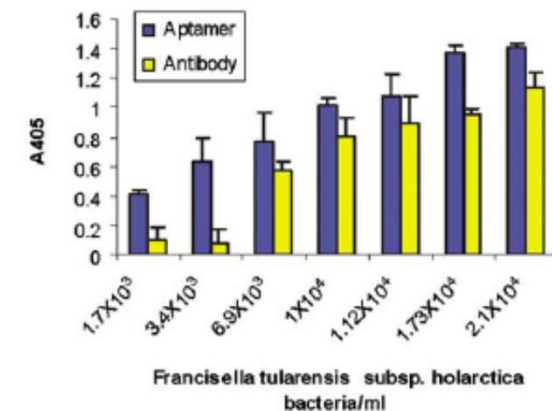


Figure 2 Tularemia bacterial antigen binding to anti-tularemia aptamer cocktail and anti-tularemia polyclonal antibodies as assessed by ALISA and ELISA using HRP activity. The assays were performed as described in 'Materials and methods'. The bacterial antigen used in the binding assay was prepared from *F. tularensis* subspecies *holarctica* (live vaccine strain). The data are plotted as OD at 405 nm vs number of bacteria/ml. Averages of triplicate measurements are shown in the figure.



SELEX: Selection of Aptamers

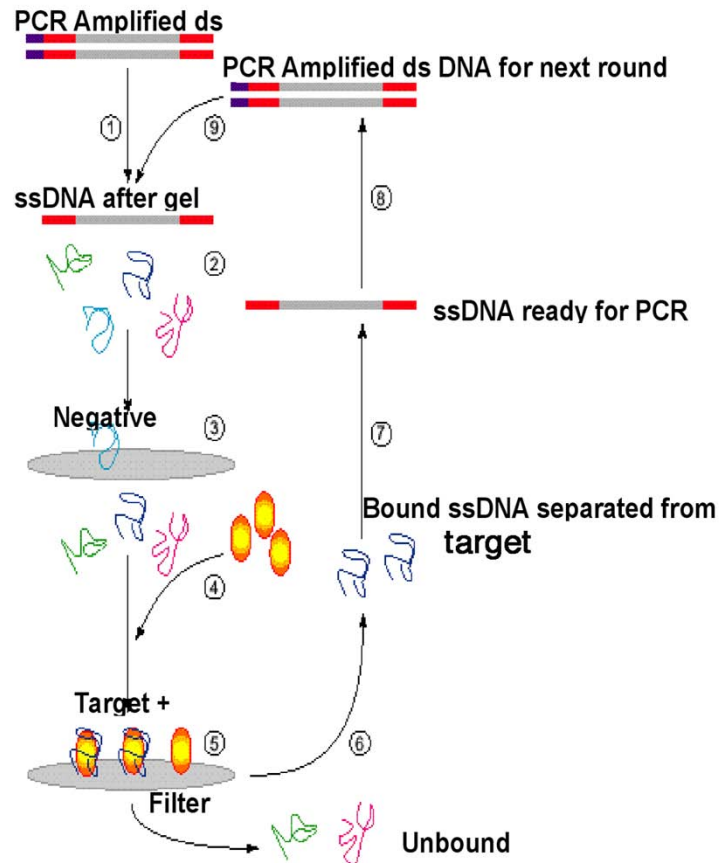

Selection begins with a library of $\sim 10^{15}$ single strands of DNA. The target is bound to a filter, and a portion of the library binds to the target. The bound strands (+) are eluted from the target by heat and amplified using PCR, with the primer for the negative strand containing biotin at its 5' end. After amplification, the DNA is denatured and the (-) strands separated from the (+) strands by passing the DNA over a streptavidin column, which retains the biotin containing (-) strands. Another round of selection is begun.

Quantum Dots

- Quantum Dots
 - Very bright
 - Resistant to photo-bleaching
 - One excitation wave length
 - Vendors
 - Evident Technologies: T1(block co-polymer) and T2 (lipid)
 - Invitrogen (polymer)

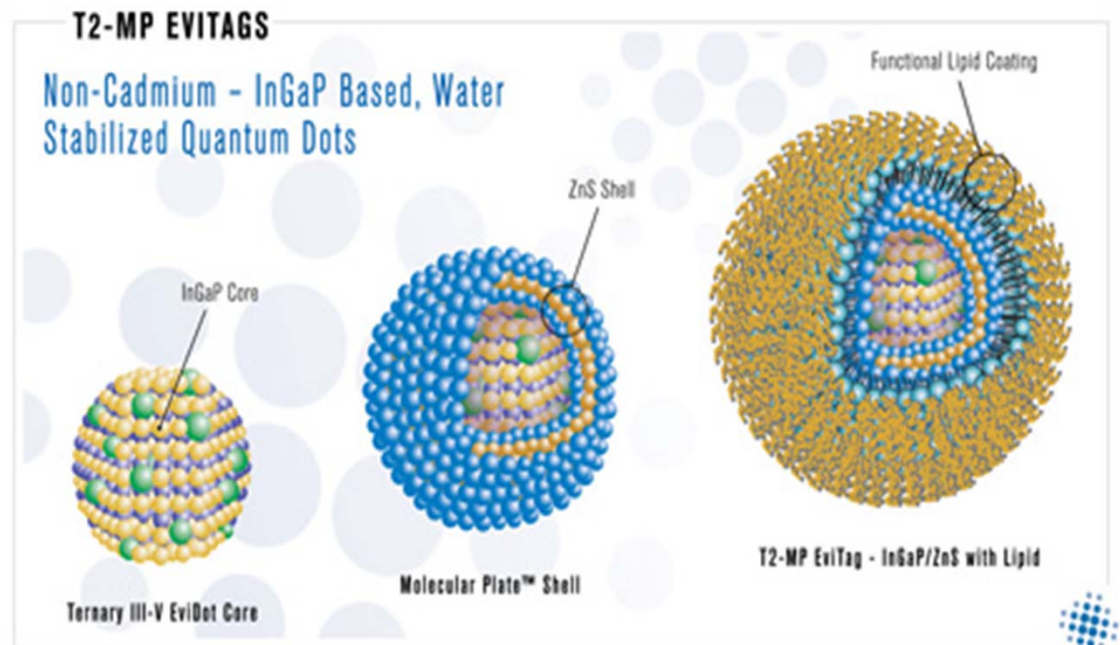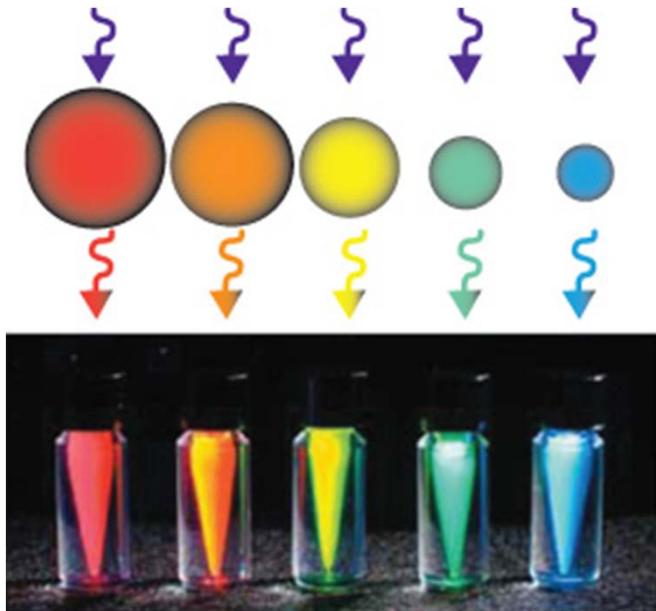


Evident Technology Web Site

Approved for public release; distribution unlimited


Immobilization of Aptamer/Qdot

Synthesis of reverse system

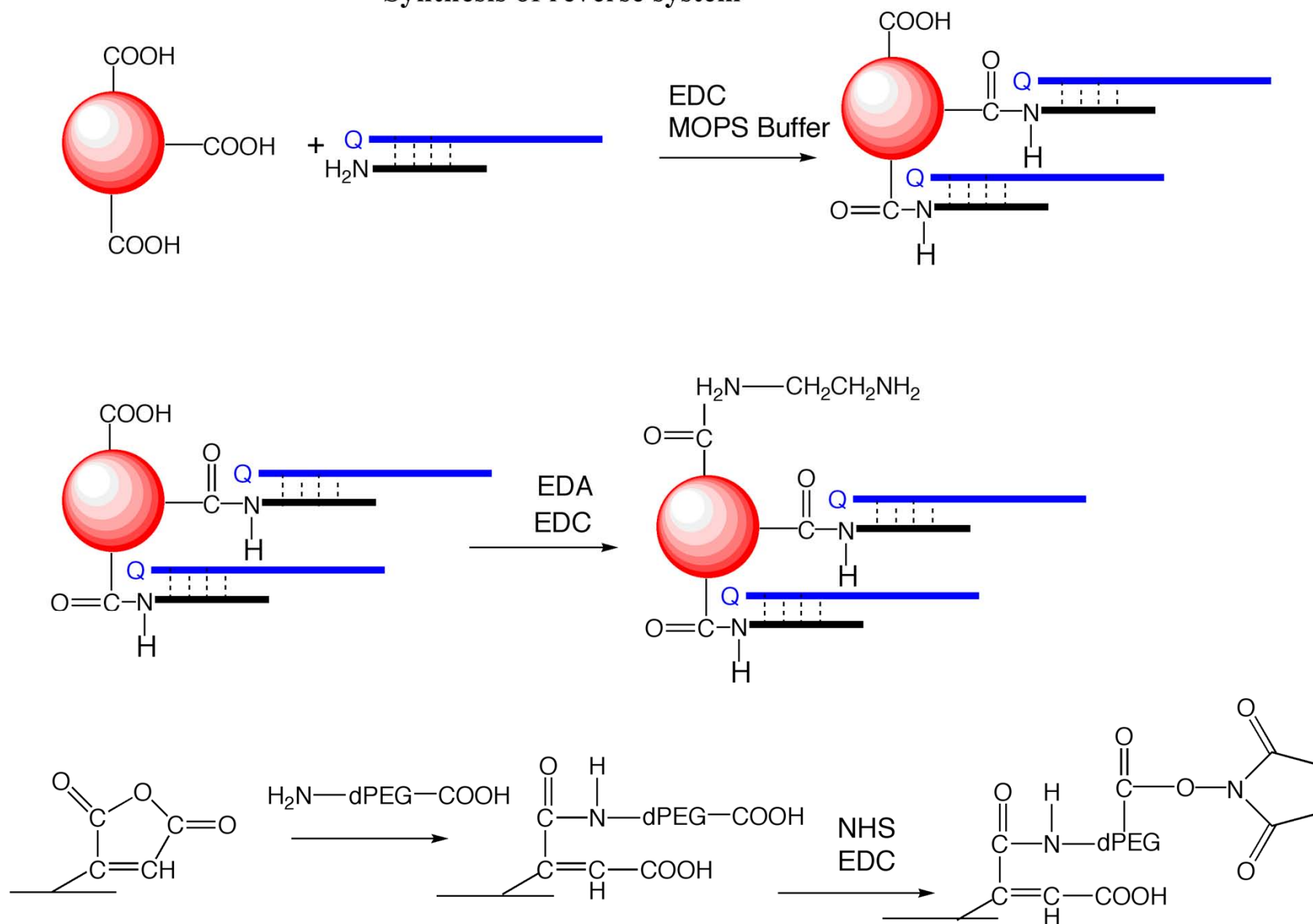


Approved for public release; distribution unlimited

Immobilization of Aptamer/Qdot

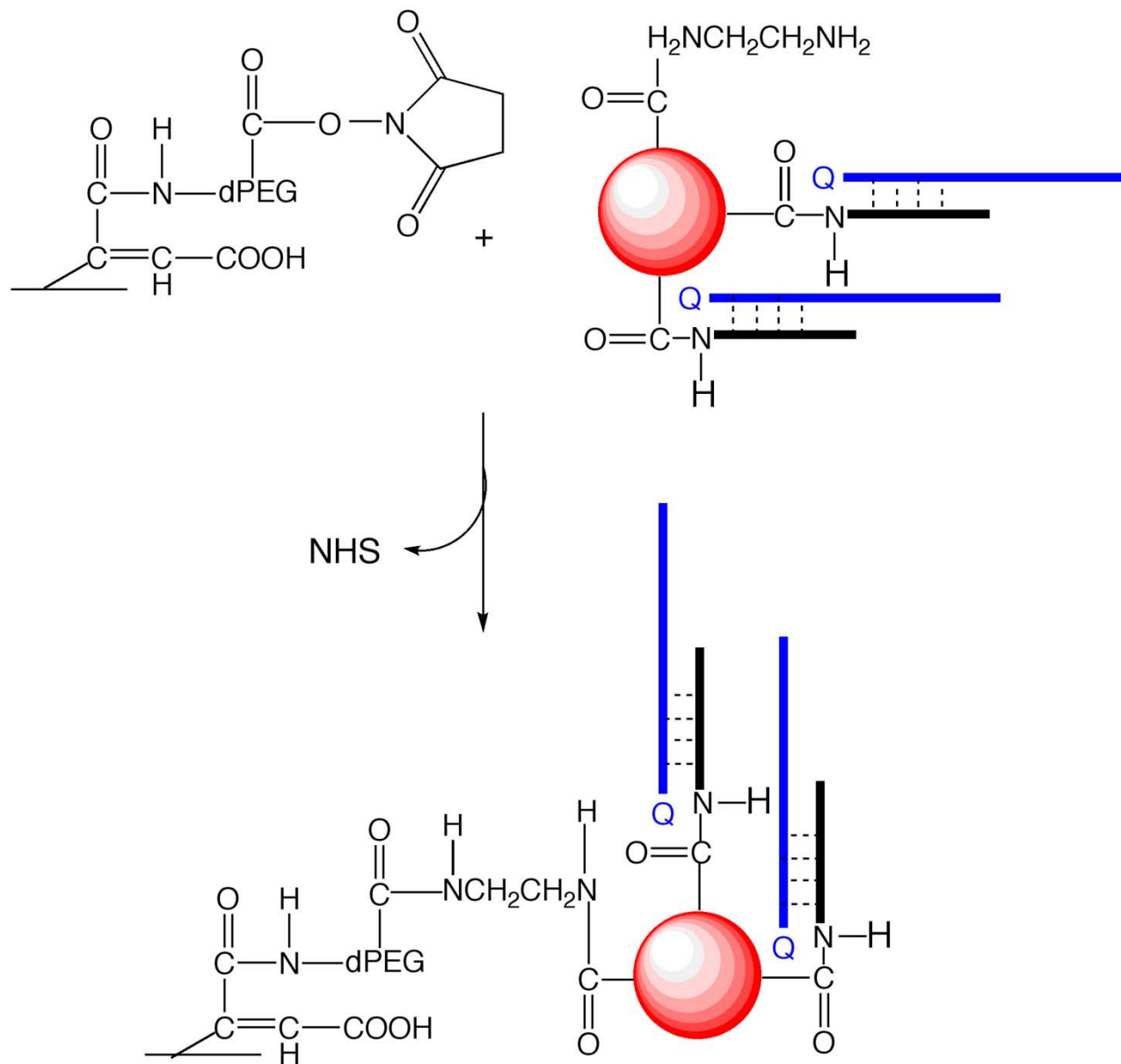

Approved for public release; distribution unlimited

Quenching/Dequenching

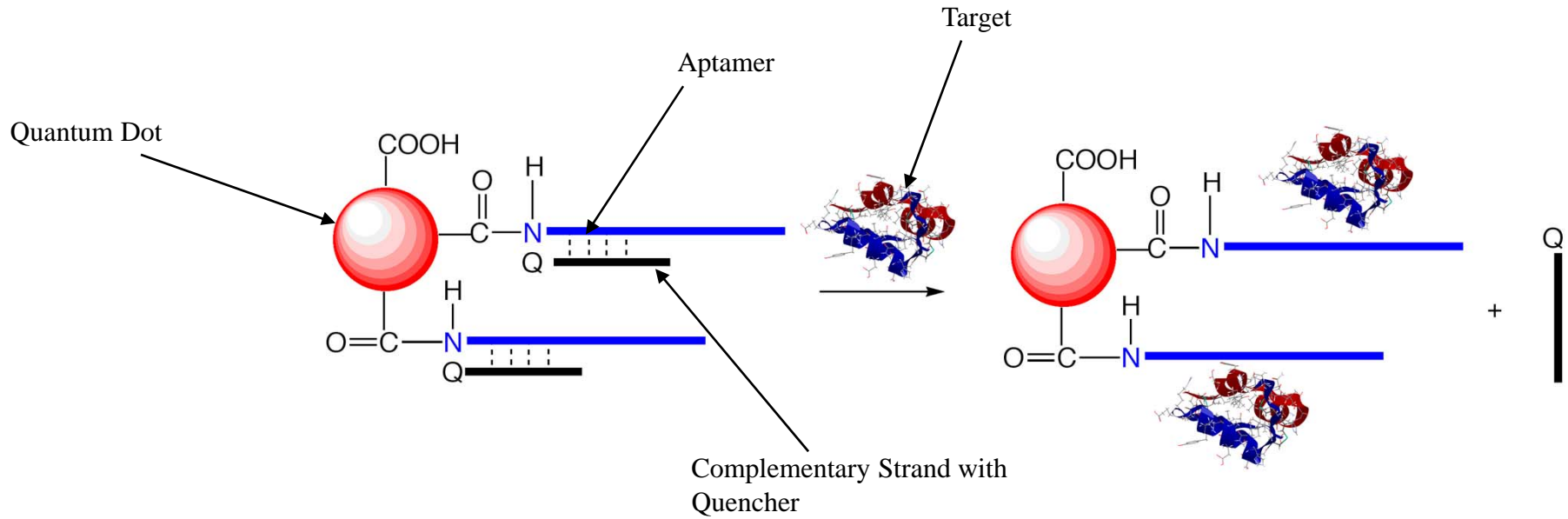


Quenched

Dequenched

Approved for public release; distribution unlimited

Possible type of Assays

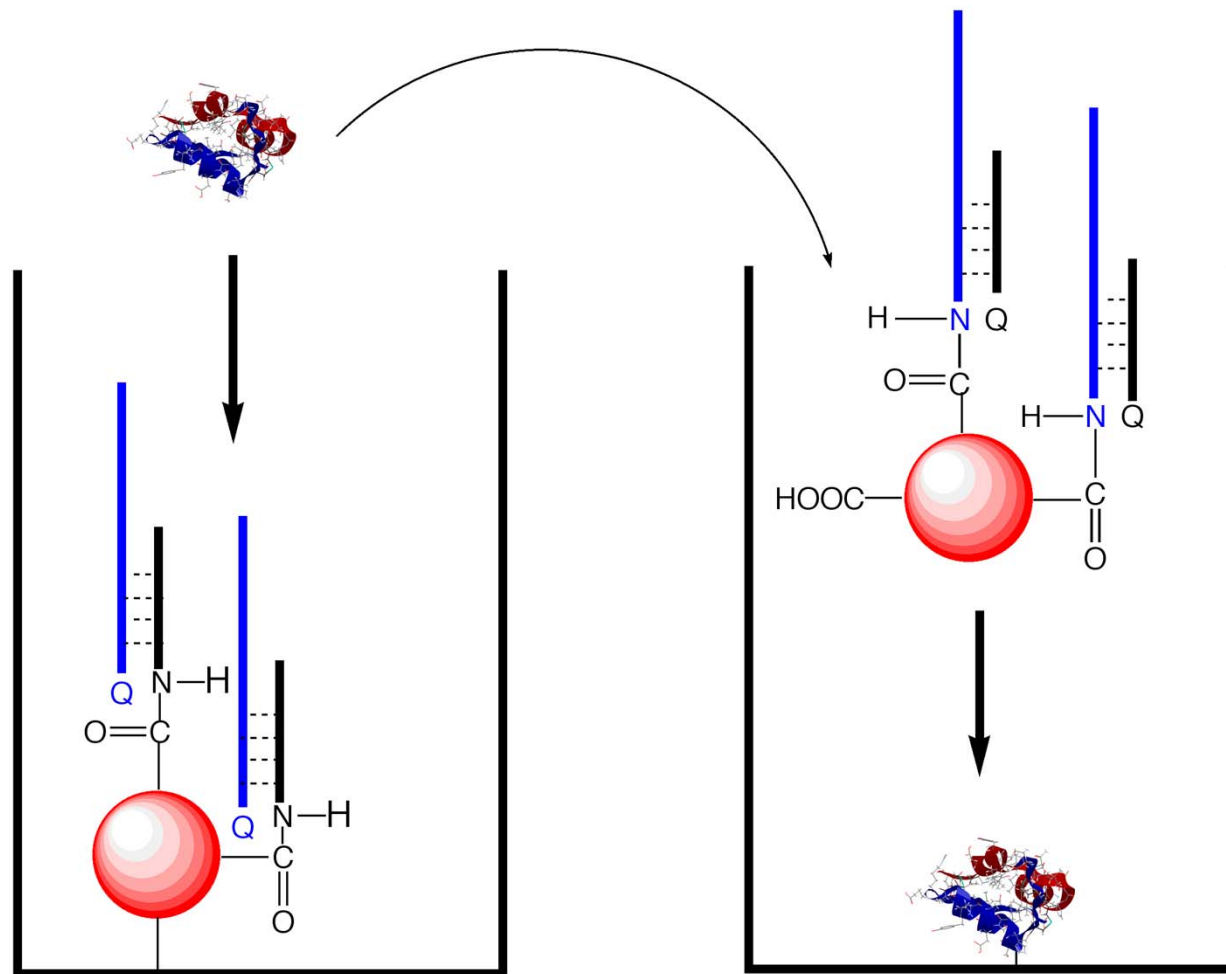


ELISA-like Assay

Competitive type Assay

Approved for public release; distribution unlimited

Aptamer/Quenched Quantum Dots Response to Shiga Toxin or Ovalbumin using SELEX DNA Aptamers


Kiel, J. L., Holwitt, E. A., and Sorola, V. K. Select Agent Recovery and Identification Using Aptamer-Linked Immobilized Sorbent Assay. Proceedings of the CB Medical Treatment Symposium, Spietz Laboratory, Switzerland **CBMTS VII** (electronic publication), 33, 1-7 (2008)

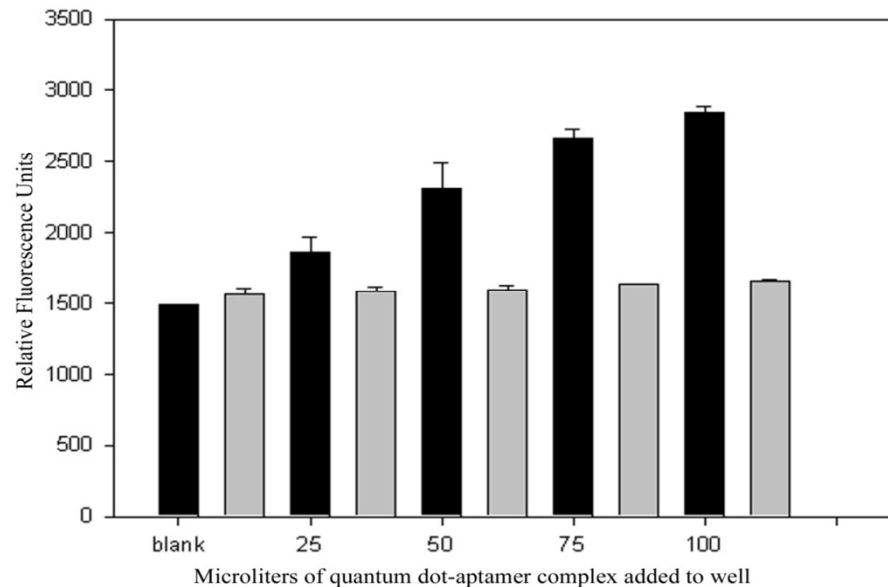

Shiga toxin: black bars
Ovalbumin: gray bars

Approved for public release; distribution unlimited

ASExpP

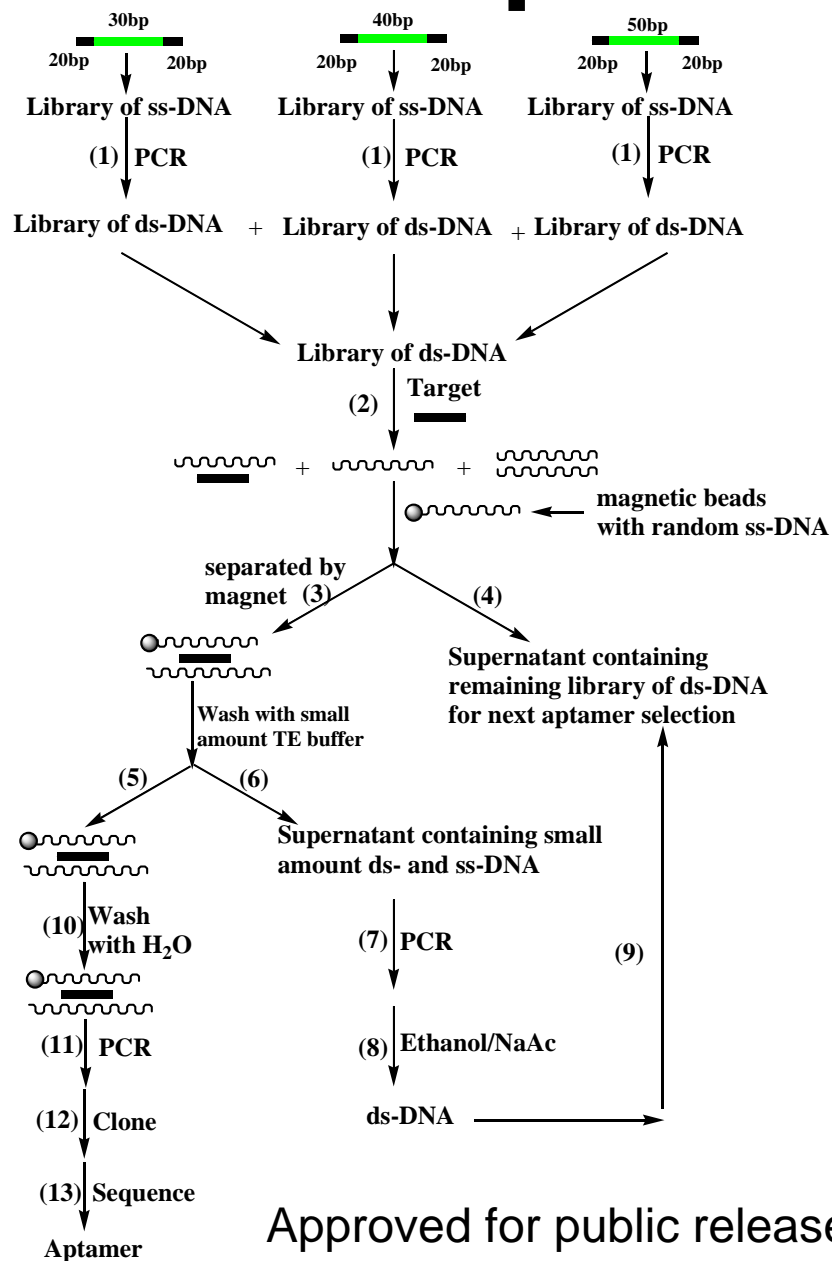

Approved for public release; distribution unlimited

Double-Strand DNA Response De-quenching Using ASExpP vs. SELEX Aptamers against Bot Tox A

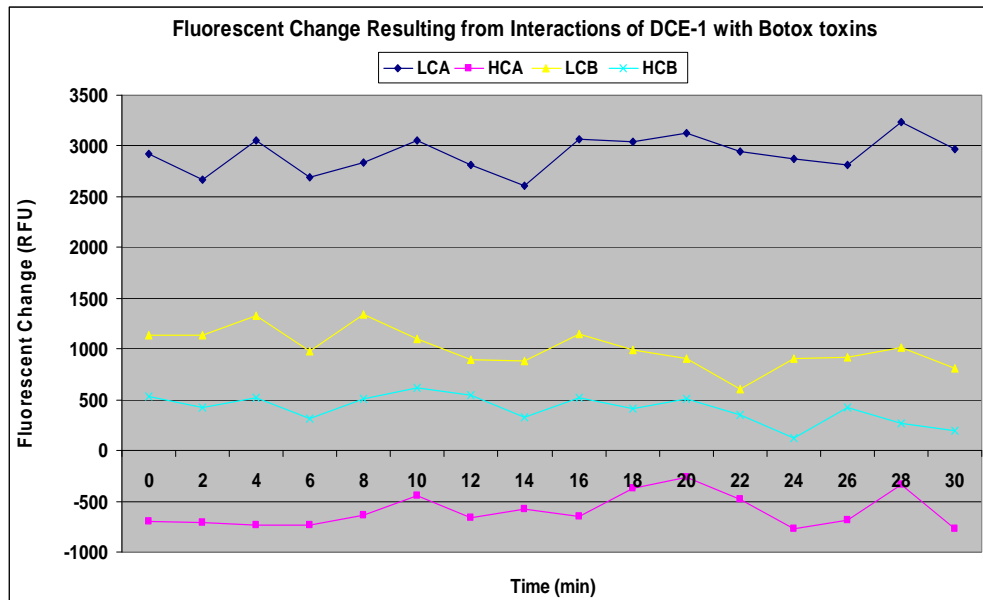


Figure A. Fluorescence change resulting from the interactions of DCE-1 (made from aptamer against BoTox, type A-light chain by **ASExpP process**) with different types of BoTox.

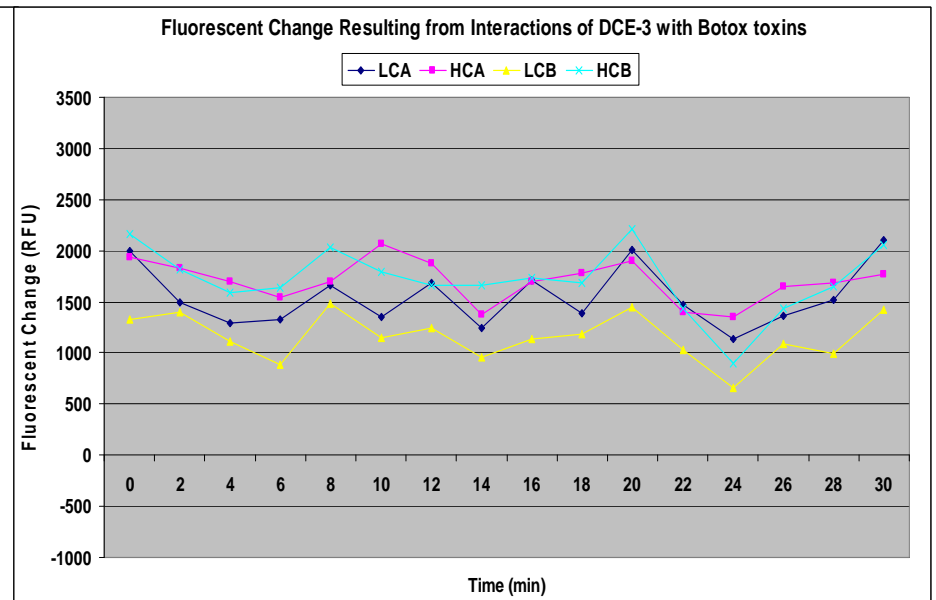


Figure B. Fluorescence change resulting from the interactions of DCE-3 (made from aptamer against BoTox, type A-light chain by **SELEX process**) with different types of BoTox.

•Fan,M., McBurnett, S. R., Andrews, C. J., Allman, A. M., Bruno, J. G., and Kiel, J. L..
Aptamer Selection Express: A Novel Method for Rapid Single-Step Selection and Sensing of
Aptamers. J Biomol Tech **19**(5), 311–319 (December 2008).

Approved for public release; distribution unlimited


Bot Tox Aptamers: **SELEX** and **ASExpP**

Selected by ASExpP against BoTox, type A-light chain (for DCE-1)
1(+). AgTCTAgAgggCCCCAgAATACACCCgACAACtAgAT ACCCATCAAAAgtCCAgCAAaggATgCaggggT
1(-). ACCCCTgCATCCTTTgCTggACTTTTgATgggTATCTA gTTgTCgggTgTATTCTggggCCCTCTAgACT
Selected by ASExpP against BoTox, type B-light chain (for DCE-2)
2(+). AgTCTAgAgggCCCCAgAATTATCCACTAgCgggAAgT AgTACATCTCACCCAgCAAaggATgCaggggT
2(-). ACCCCTgCATCCTTTgCTgggTgAgATgTACTACTTCC CgCTAgTggATAATTCTggggCCCTCTAgACT
Selected by SELEX against BoTox, type A-light chain (for DCE-3)
3(+). CATCCgTCACACCTgCTCTggggATgTgTggTgTTggCT CCCgTATCAAgggCgAATTCT
3(-). gTAggCagTgTggACgAgACCCCTACACACCACAACC gAgggCATAgTTCCCgCTTAAGa
Selected by SELEX against BoTox Holotoxin (for DCE -4)
4(+). CATCCgTCACACCTgCTCTgCTATCACATgCCTgCTg AAgTggTgTTggCTCCCgTATCA
4(-). gTAggCagTgTggACgAgACgATAgTgTACggACgACTTC ACCACAACCgAgggCATAgT

Approved for public release; distribution unlimited

Magnetic Nanoparticle or Microparticle and Quantum Dot Separation

Internal Standard Subtractive Ratio Assay Method

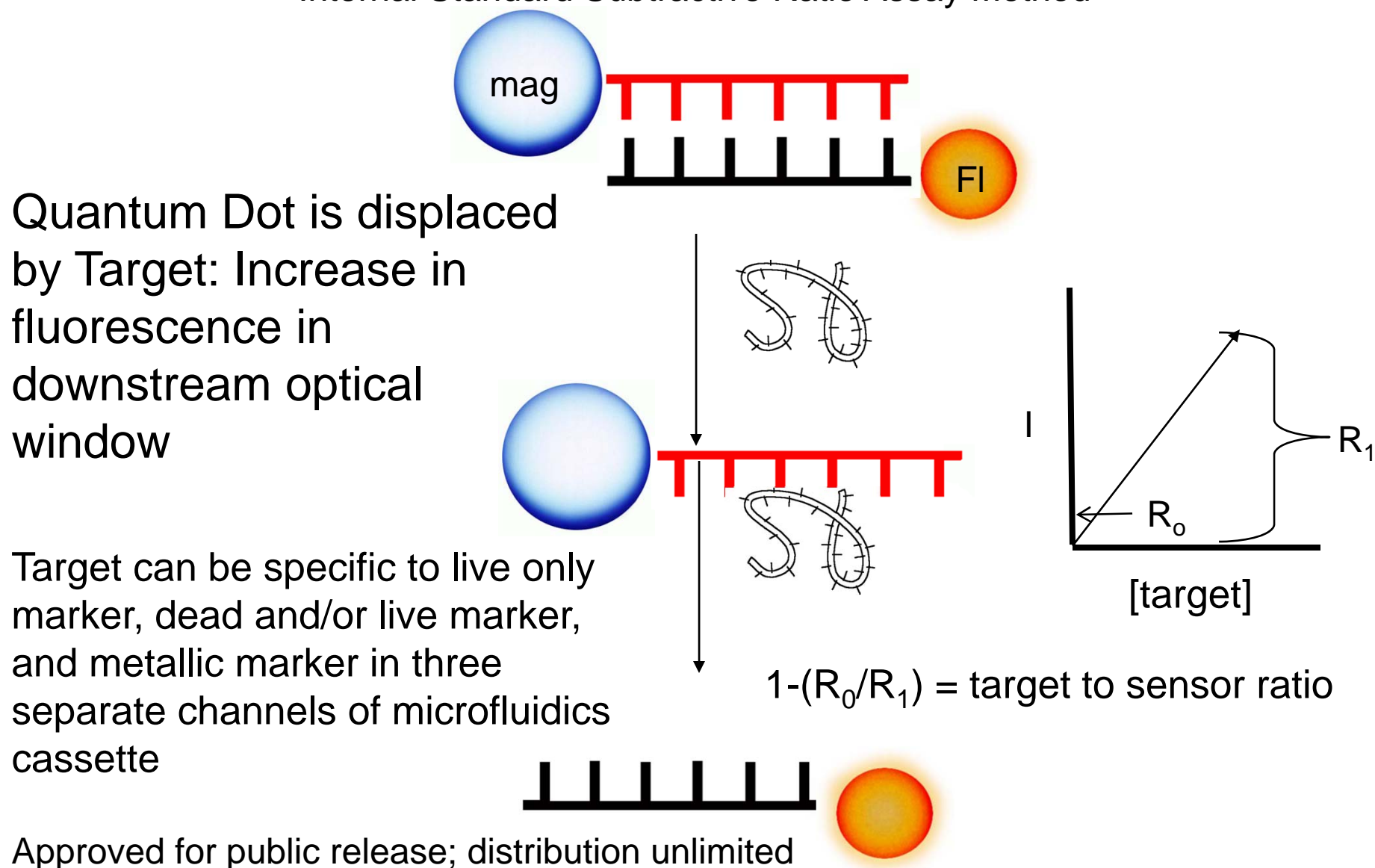

Aptamer Based Agent Detection

GOAL: Man portable detection of biological agents in the field

The Portable Test Laboratory has been flown on the International Space Station



Identifies positive sample in field allowing for further analysis in a controlled location !




The Portable Test Laboratory has been tested in conditions of extreme heat and cold


Approved for public release; distribution unlimited

AFRL/CRL Device


Currently Marketed Portable Test Lab

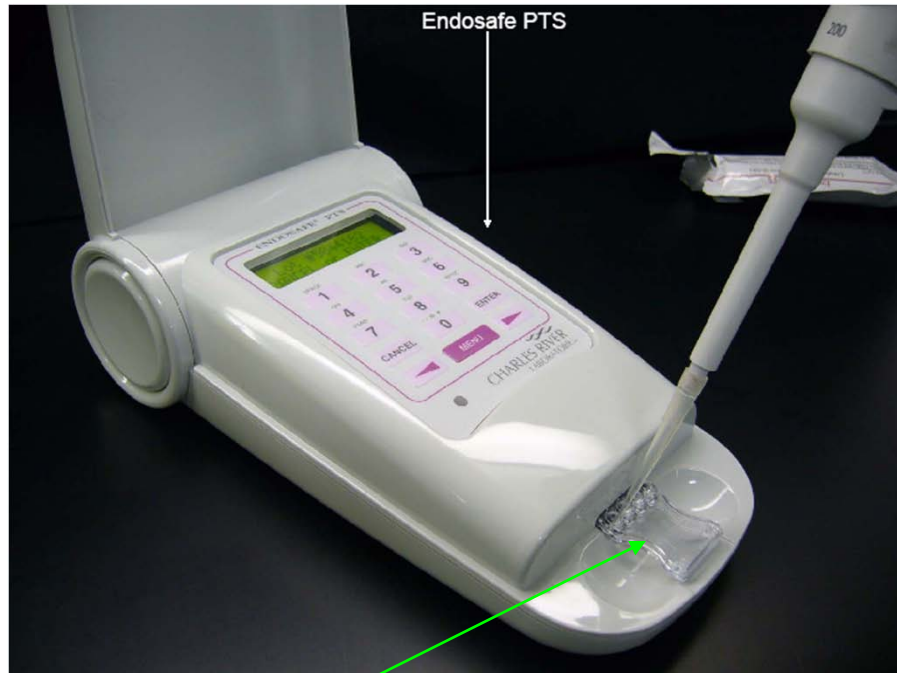

Cassette loaded into
Portable Test
Laboratory

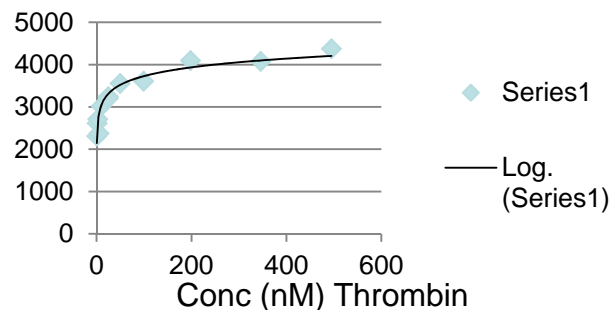


- Cartridge provides ability to retrieve samples for identification and analysis of substance
- Device tests and identifies the contained specimen
- Current system measures fluorescence
- 9.25x4.625x2.50 inches
- Battery Operation for 4 hrs
- ~ 2 lbs

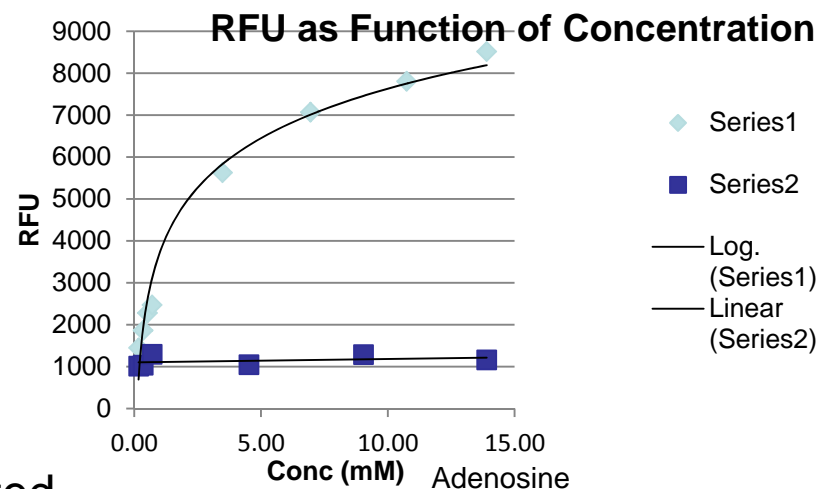


Approved for public release; distribution unlimited

GE GRC Concepts: Complementary Sensing Structures


Bare electrode structures

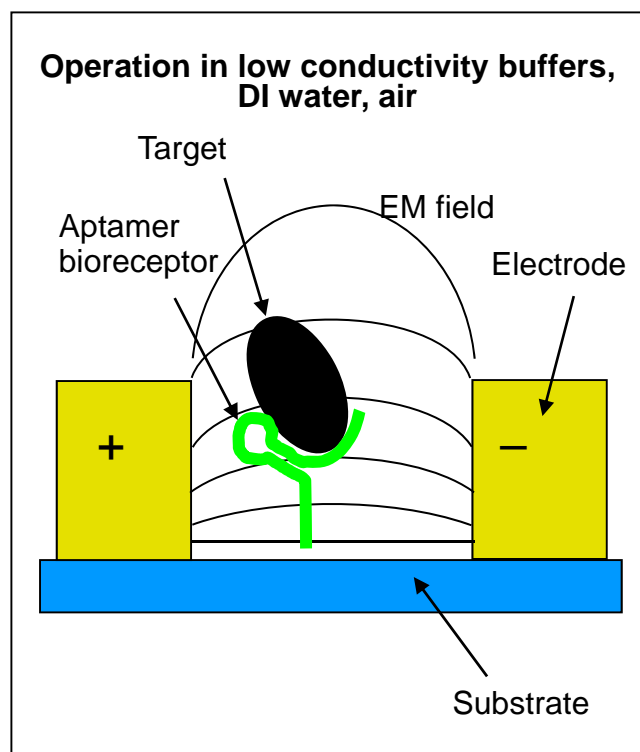


SiO₂-coated electrode structures

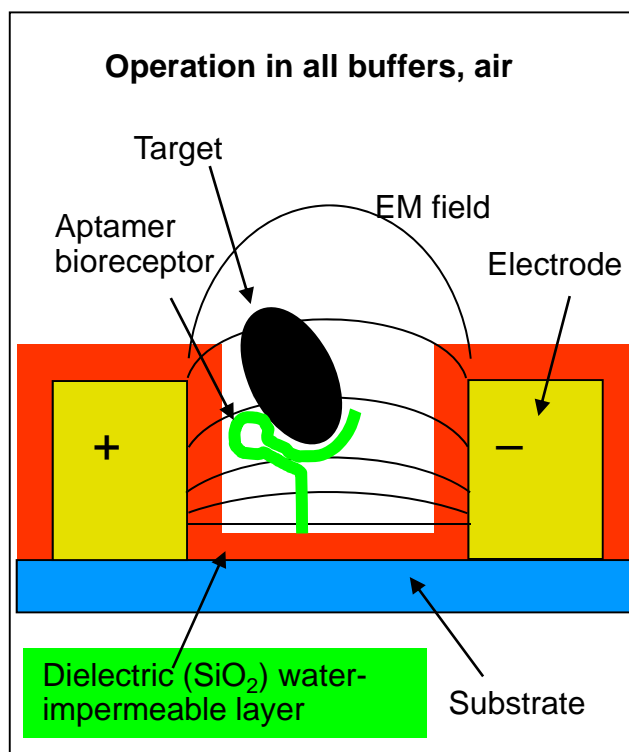

3D model of EM field

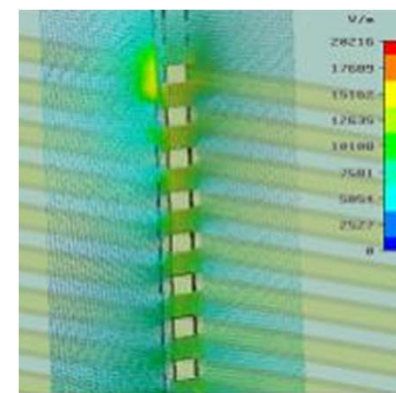

Nanofabricated sensor

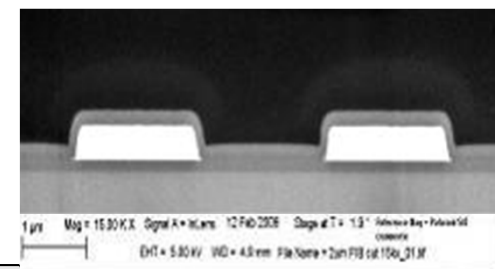


Detection of changes in capacitance and resistance of sensing gap between electrodes provides improved sensitivity and stability and rejects interferences

R. A. Potyrailo, C. Surman, R. Chen, S. Go, K. Dovidenko, and W. G. Morris, General Electric Company, Global Research Center, Niskayuna, NY, USA; E. Holwitt, V. Sorola, and J. L. Kiel, Air Force Research Lab RHPC, Brooks City-Base, TX, USA. Label-Free Biosensing Using Passive Radio-Frequency Identification (RFID) Sensors. 15th International Conference on Solid-State Sensors, Actuators and Microsystems - Transducers'09, 21-25, June 2009, Denver, Colorado USA

Approved for public release; distribution unlimited

GE GRC: Two Designs of RFID Sensing Electrodes

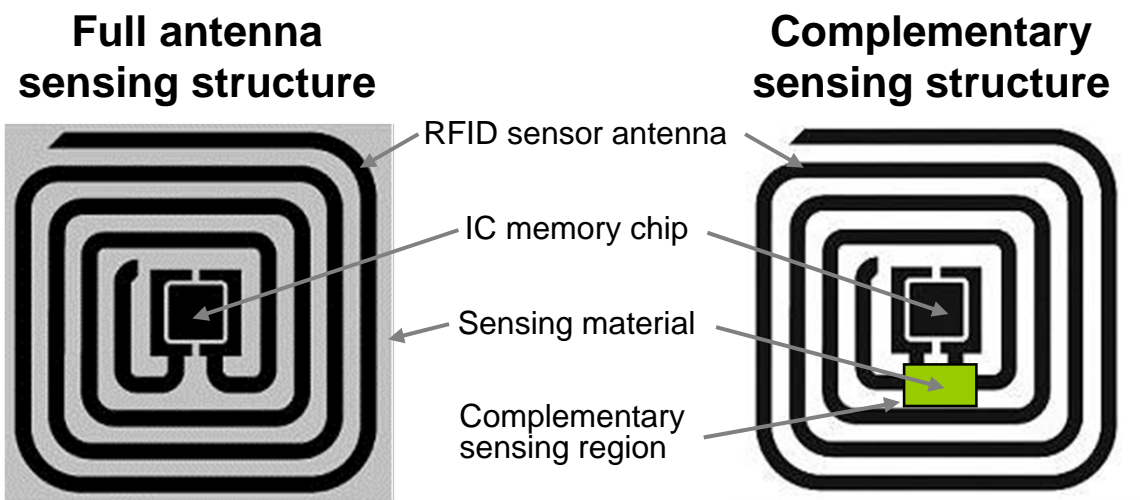

Pros

- Simple design
- Ease of fabrication

Cons

- Reagent cost

Pros

- Smaller sensing area
- Ease to deposit sensing films
- Highest sensitivity

Cons


- Medium fabrication difficulty

R. A. Potyrailo, C. Surman, R. Chen, S. Go, K. Dovidenko, and W. G. Morris, General Electric Company, Global Research Center, Niskayuna, NY, USA; E. Holwitt, V. Sorola, and J. L. Kiel, Air Force Research Lab RHPC, Brooks City-Base, TX, USA. Label-Free Biosensing Using Passive Radio-Frequency Identification (RFID) Sensors. 15th International Conference on Solid-State Sensors, Actuators and Microsystems - Transducers'09, 21-25, June 2009, Denver, Colorado USA


Approved for public release; distribution unlimited

GE GRC: Analysis of K_d for Thrombin Aptamer

Reference	Aptamer use	Sensing format	K_d
Ostatna, V. et al. <i>Anal. Bioanal. Chem.</i> 2008, 391(5), 1861	3'-biotin or 3'-SH onto avidin, Au and dendrimer surfaces	SPR	40-385nM
Potyraiilo, R., et al. <i>Anal. Chem.</i> , 1998, 70, 3419.	3'-C7 glass slide immobilized	ATR-fluorescence anisotropy	47nM
Lee, M.; Walt, D. R., <i>Anal. Biochem.</i> 2000, 282, (1), 142.	5'-NH-C6 on silica beads	Competition, fluorescence	300nM
GE GRC (THIS WORK)	Functionalized aptamers on glass slide	Fluorescence	10-100nM
Li, J. J. et al. <i>Biochem. Biophys. Res. Comm.</i> 2002, 292, (1), 31.	Solution	Molecular beacon, fluorescence	5.20 ± 0.49 nM
Hamaguchi, N., et al. <i>Anal. Biochem.</i> 2001, 294, (2), 126.	Solution	Molecular beacon, fluorescence	10nM
Tasset, D. M. et al. <i>J. Mol. Biol.</i> 1997, 272, (5), 688.	Solution	Filter binding	75-100nM

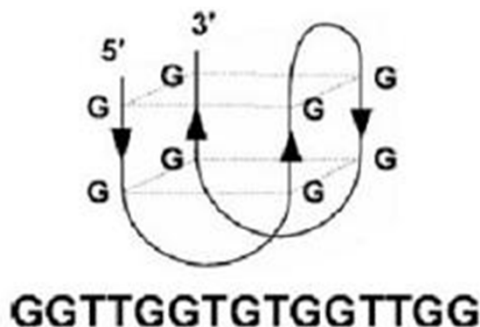



Thrombin aptamer: Bock et al., *Nature* 355: 564-566, 1992.

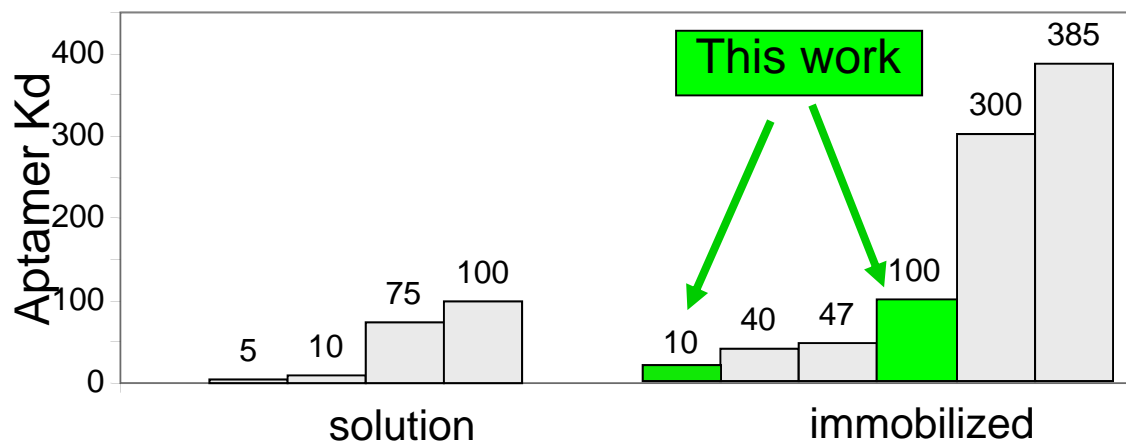

R. A. Potyraiilo, C. Surman, R. Chen, S. Go, K. Dovidenko, and W. G. Morris, General Electric Company, Global Research Center, Niskayuna, NY, USA; E. Holwitt, V. Sorola, and J. L. Kiel, Air Force Research Lab RHPC, Brooks City-Base, TX, USA. Label-Free Biosensing Using Passive Radio-Frequency Identification (RFID) Sensors. 15th International Conference on Solid-State Sensors, Actuators and Microsystems - Transducers'09, 21-25, June 2009, Denver, Colorado USA

Approved for public release; distribution unlimited

GE GRC: *Spore detection:* Characterization of 2-D Bare Nanogap FIB-fabricated Electrodes

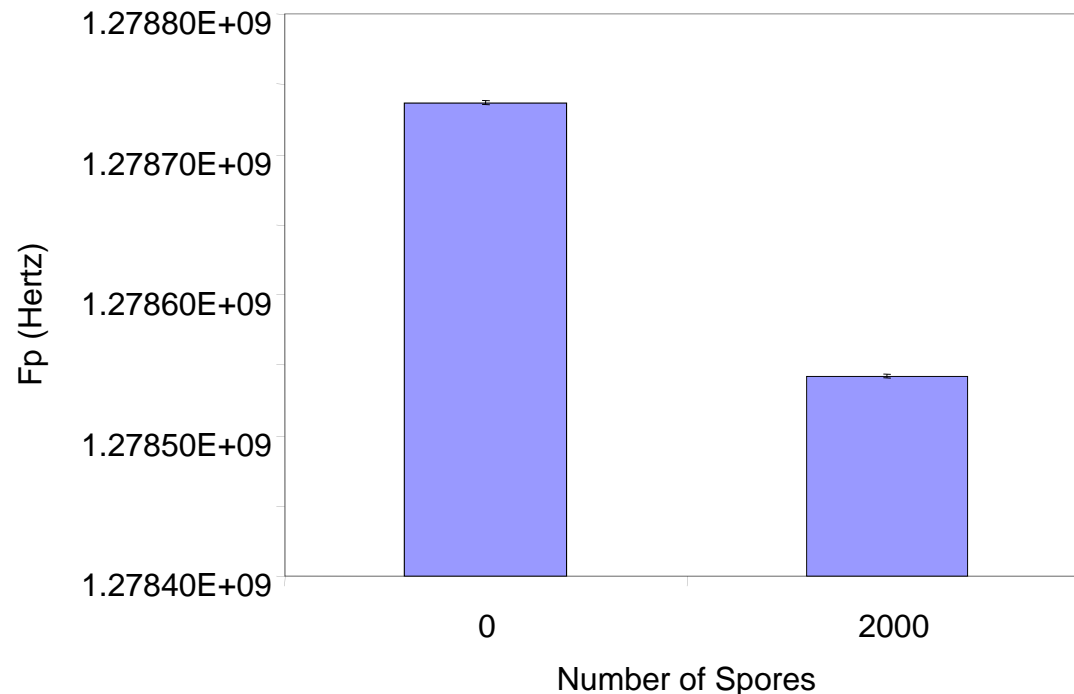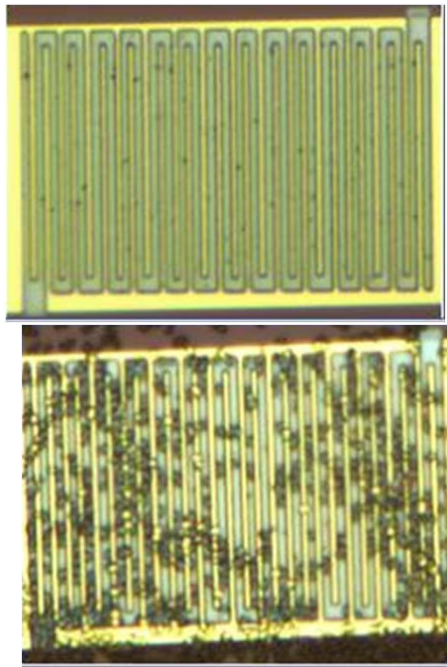

Detection limit of BG spores = 35 spores
Most techniques except for culture (1 spore)
detect a minimum of 100-100,000 spores

Emerging Exotic Pathogens: Heartwater and Viper Plague



Amblyomma species responsible for transmission of Ehrlichia ruminantium (Photo courtesy of APHIS-USDA)


- **Heartwater**
 - Tick-borne disease: *Amblyomma variegatum*, *A. hebraeum*, *A. lepidum*, *A. maculatum*, other Amblyomma tick carriers
 - Causal agent: *Cowdria ruminantium*, now *Ehrlichia ruminantium*
- **Imminent threat to Western Hemisphere**
 - Mortality in cattle and other ruminants: excess of 70%
 - Has been found in African spurred tortoises (*Geochelone sulcata*) and leopard tortoises (*Geochelone pardalis*)
 - Is now in Caribbean Islands
 - Antigua
 - Guadeloupe
 - Marie Galante
 - Perhaps Cuba
- **Viper Plague, a mimic of heartwater, and associated ticks entered the USA in 2002**
- **VP rickettsia was isolated in viper cells and propagated in turtle cells, but also infects bovine endothelial cells, and human cells (HeLa)**

Kiel, J. L. , Alarcon, R. M., Parker, J. L., Vivekananda, J., Gonzalez, Y. B., Stribling, L. J. V., and Andrews, C., Emerging Tick-Borne Disease in African Vipers Caused by a Cowdria-Like Organism, Ann. N.Y. Acad. Sci. 1081: 434-442, 2006

Approved for public release; distribution unlimited

Molecular Biology Confusion (Standard Diagnostic PCR) Between Heartwater and Viper Plague

Kiel, J.L., Gonzalez, Y., Parker, J.E.,
Andrews, C., Martinez, D., Vacheiry, N.,
LeFrancois, T. Viral association with the
elusive rickettsia of viper plague from
Ghana, West Africa. *Annals of the New
York Academy of Sciences* 1149, 318-321
(2008).

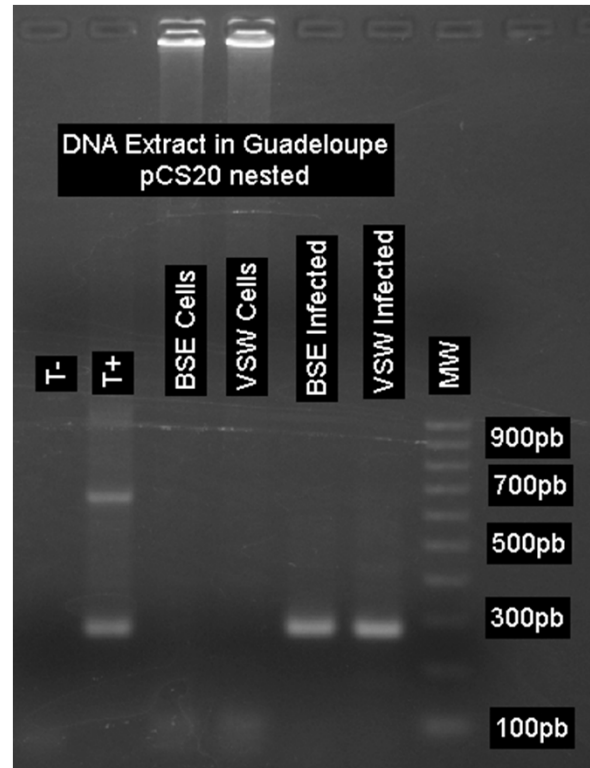

**Nested PCR pCS20:
AB128/129/130**

PCR products sent for sequencing:

PCR pCS20F-HpCS20R: 750pb instead of 1100pb

Nested pCS20: 280pb like *Ehrlichia ruminantium* (heartwater agent)

VSW and BSE ER: Viper spleen and Bovine endothelial
cells infected


VSWC & BSEC: uninfected cells

T+= positive control DNA from *Ehrlichia ruminantium*

MW= ladder 100pb


Approved for public release; distribution unlimited

Centrifuged Bovine Endothelial Cell Supernatant Showing Rickettsia (requires many large culture flasks to accumulate this number of rickettsia)

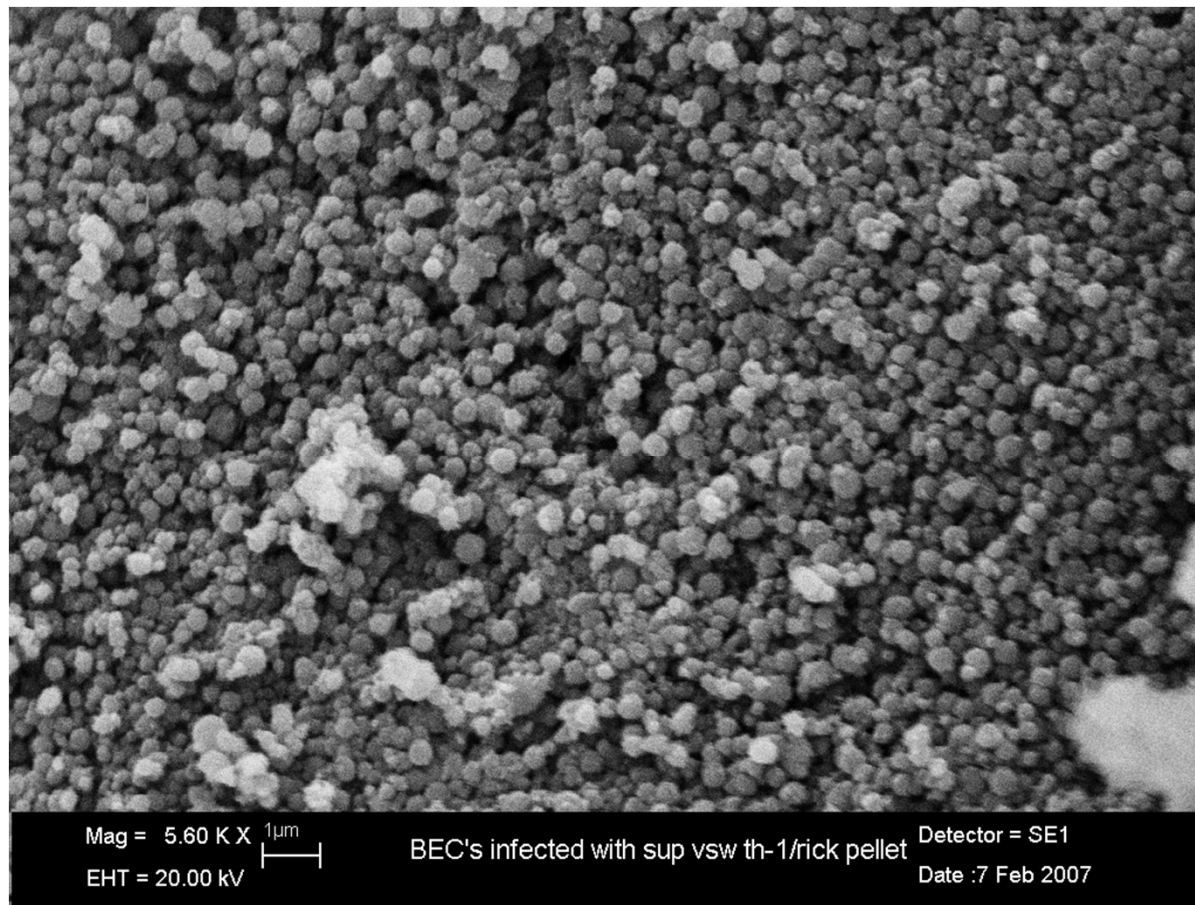


Approved for public release; distribution unlimited

Bovine Endothelial Cells: Infected with VP Showing a Hidden Type D Immunosuppressive Retrovirus associated with the Disease Compared to Human Type D Retrovirus: Could VP Type D Retrovirus be the First Biological Vector (Tick) Transmitted Retrovirus?

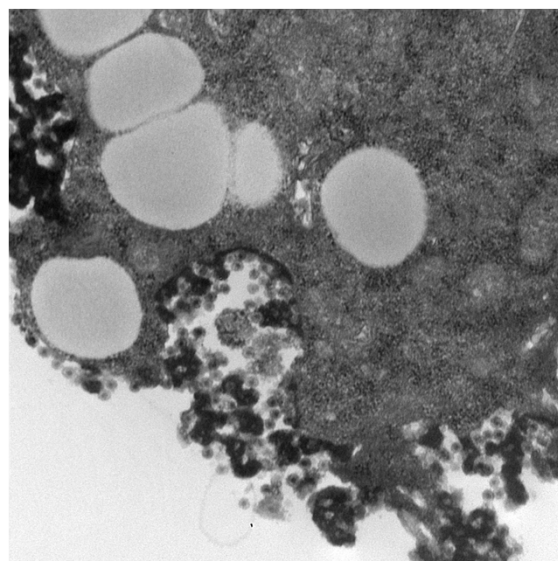

BROOKS CITY BASE.020.tif
BOVINE
Print Mag: 16500x @ 8. in
10:02 09/20/06
Microscopist: EM


500 nm

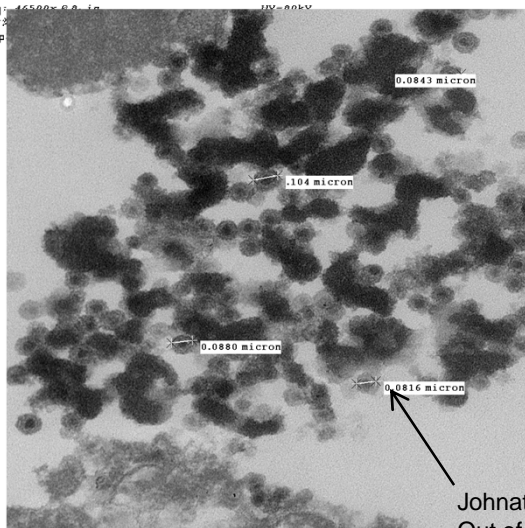

BROOKS CITY BASE.011.tif
BOVINE
Print Mag: 93100x @ 8. in
9:46 09/20/06
Microscopist: EM

500 nm

HV=80kV
Direct Mag: 60000x
UTHSC-SA Pathology

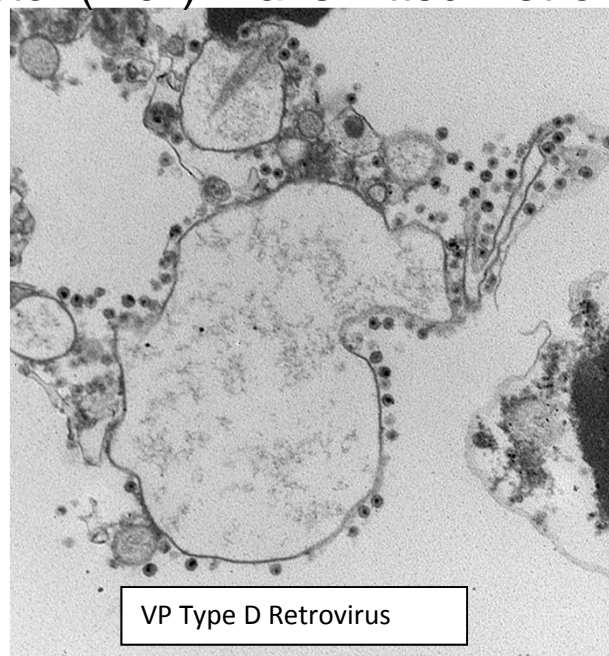


VP Type D Retrovirus

BROOKS CITY BASE.196.tif
BULL SERUM
Print Mag: 38800x @ 8. in
10:24 01/23/07
Microscopist: EM

500 nm

HV=80kV
Direct Mag: 25000x
UTHSC-SA Pathology

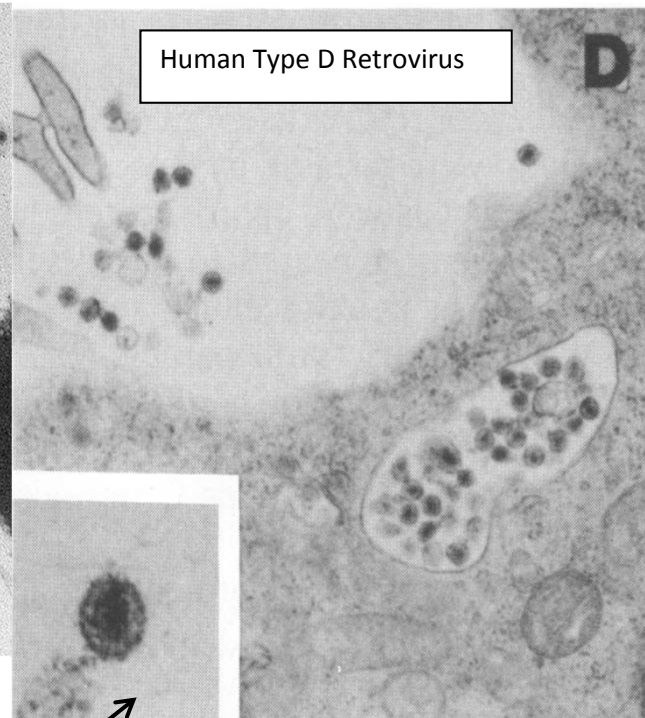

Human Type D Retrovirus

Isolation of a Type D Retrovirus from B-Cell Lymphomas of a Patient with AIDS

ROBERT C. BOHANNON,¹² LAWRENCE A. DONEHOWER,¹ AND RICHARD J. FORD^{2*}

Division of Molecular Virology, Baylor College of Medicine,¹ and Department of Molecular Pathology, University of Texas M. D. Anderson Cancer Center,² Houston, Texas 77030


Received 29 November 1990/Accepted 23 July 1991

JOURNAL OF VIROLOGY, Nov. 1991, p. 5663-5672

0022-538X/91/115663-10\$02.00/0



Copyright © 1991, American Society for Microbiology

Johnathan L. Kiel, Yvette Gonzalez, Ishmael I. Rosas and David F. Vela,
Out of Africa: Do Viruses Play a Role in the Emergence of New Rickettsial Diseases?
Presentation at the 5th International Meeting on Rickettsiae and Rickettsial Diseases. Marseille, France. May 2008

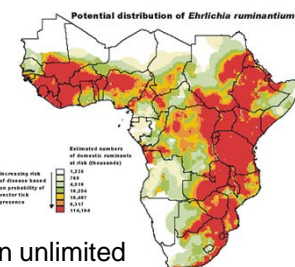

Approved for public release; distribution unlimited

New Retrovirus Infects a Wider Host Range than VSW Virus

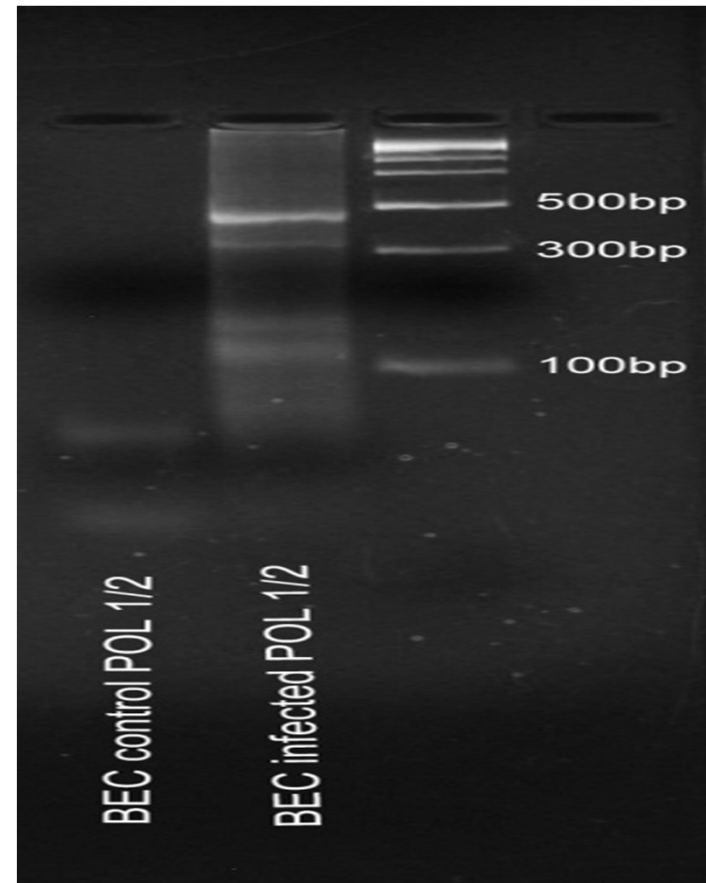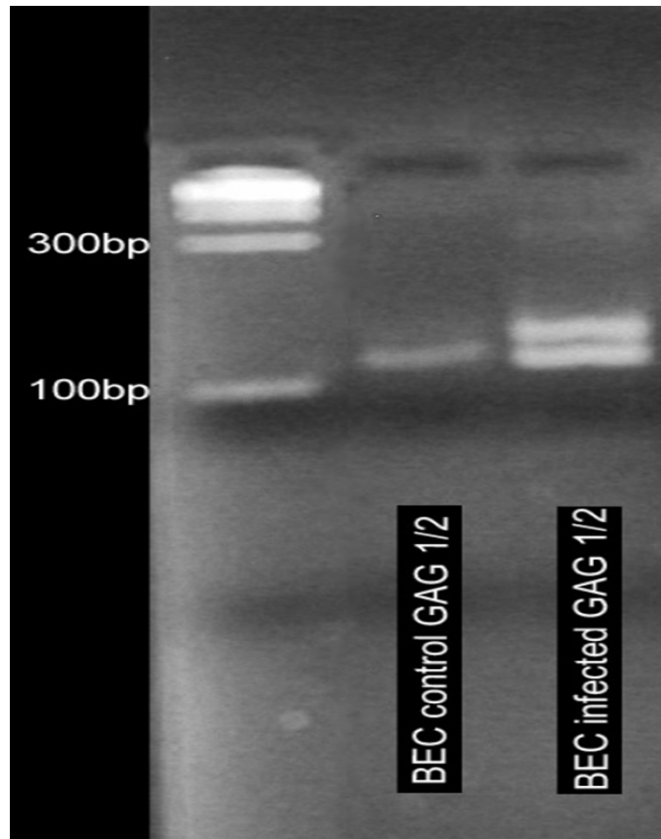

J. Kiel, Y. Gonzalez, J. Parker, C. Andrews, D. Martinez, N. Vachieri, T. Lefrancois, ANYAS 1149: 318-321 (2008)

Approved for public release; distribution unlimited

It started in the Garden of Eden: The Serpent and the Tree of the Knowledge of Good and Evil

Snakes as agents of evolutionary change in primate brains

Journal of Human Evolution 51 (2006) 1-35

Lynne A. Isbell*

Department of Anthropology, University of California, Davis, CA 95616, USA

Received 6 March 2004; accepted 28 December 2005

Abstract

Current hypotheses that use visually guided reaching and grasping to explain orbital convergence, visual specialization, and brain expansion in primates are open to question now that neurological evidence reveals no correlation between orbital convergence and the visual pathway in the brain that is associated with reaching and grasping. **An alternative hypothesis proposed here posits that snakes were ultimately responsible for these defining primate characteristics. Snakes have a long, shared evolutionary existence with crown-group placental mammals and were likely to have been their first predators.** Mammals are conservative in the structures of the brain that are involved in vigilance, fear, and learning and memory associated with fearful stimuli, e.g., predators. Some of these areas have expanded in primates and are more strongly connected to visual systems. However, primates vary in the extent of brain expansion. This variation is coincident with variation in evolutionary co-existence with the more recently evolved venomous snakes. Malagasy prosimians have never co-existed with venomous snakes, New World monkeys (platyrrhines) have had interrupted co-existence with venomous snakes, and Old World monkeys and apes (catarrhines) have had continuous coexistence with venomous snakes. The koniocellular visual pathway, arising from the retina and connecting to the lateral geniculate nucleus, the superior colliculus, and the pulvinar, has expanded along with the parvocellular pathway, a visual pathway that is involved with color and object recognition. I suggest that expansion of these pathways co-occurred, with the koniocellular pathway being crucially involved (among other tasks) in pre-attentional visual detection of fearful stimuli, including snakes, and the parvocellular pathway being involved (among other tasks) in protecting the brain from increasingly greater metabolic demands to evolve the neural capacity to detect such stimuli quickly. A diet that included fruits or nectar (though not to the exclusion of arthropods), which provided sugars as a neuroprotectant, may have been a required preadaptation for the expansion of such metabolically active brains. Taxonomic differences in evolutionary exposure to venomous snakes are associated with similar taxonomic differences in rates of evolution in cytochrome oxidase genes and in the metabolic activity of cytochrome oxidase proteins in at least some visual areas in the brains of primates. Raptors that specialize in eating snakes have larger eyes and greater binocularity than more generalized raptors, and provide non-mammalian models for snakes as a selective pressure on primate visual systems. **These models, along with evidence from paleobiogeography, neuroscience, ecology, behavior, and immunology, suggest that the evolutionary arms race begun by constrictors early in mammalian evolution continued with venomous snakes. Whereas other mammals responded by evolving physiological resistance to snake venoms, anthropoids responded by enhancing their ability to detect snakes visually before the strike.**

© 2006 Elsevier Ltd. All rights reserved.

Keywords: Primate origins; Prosimians; Anthropoids; Viperids; Elapids; Koniocellular; Superior colliculus; Pulvinar; Parvocellular; Cytochrome oxidase; Pre-attention

Ancient contact between venomous snakes and catarrhines is suggested by an endogenous retrovirus in the Asian Russell's viper (*Vipera russellii*), which is more closely related to catarrhine type D retroviruses (Mason-Pfizer monkey virus and langur endogenous virus) than a platyrrhine type D retrovirus (squirrel monkey retrovirus) (Andersen et al., 1979). Endogenous viruses are often evolutionarily very old (Johnson and Coffin, 1999; van der Kuyl et al., 2000). Given our knowledge that transfer of nonhuman primate retroviruses to humans typically involves physical contact (Weiss and Wrangham, 1999; Wolfe et al., 2004), it is conceivable that the retrovirus was transferred to the Russell's viper upon contact with a catarrhine primate millions of years ago.




Approved for public release; distribution unlimited

General Spotted Fever or Typhus targeting by anti-OX-19 Antigen Aptamer-coated Particles (Dr Fan by ASEXP)



Visible and UV Light Photomicrographs after
OX-19 Fluorescent Antibody Treatment

Control VH2 Cells

RLO Infected VH2 Cells

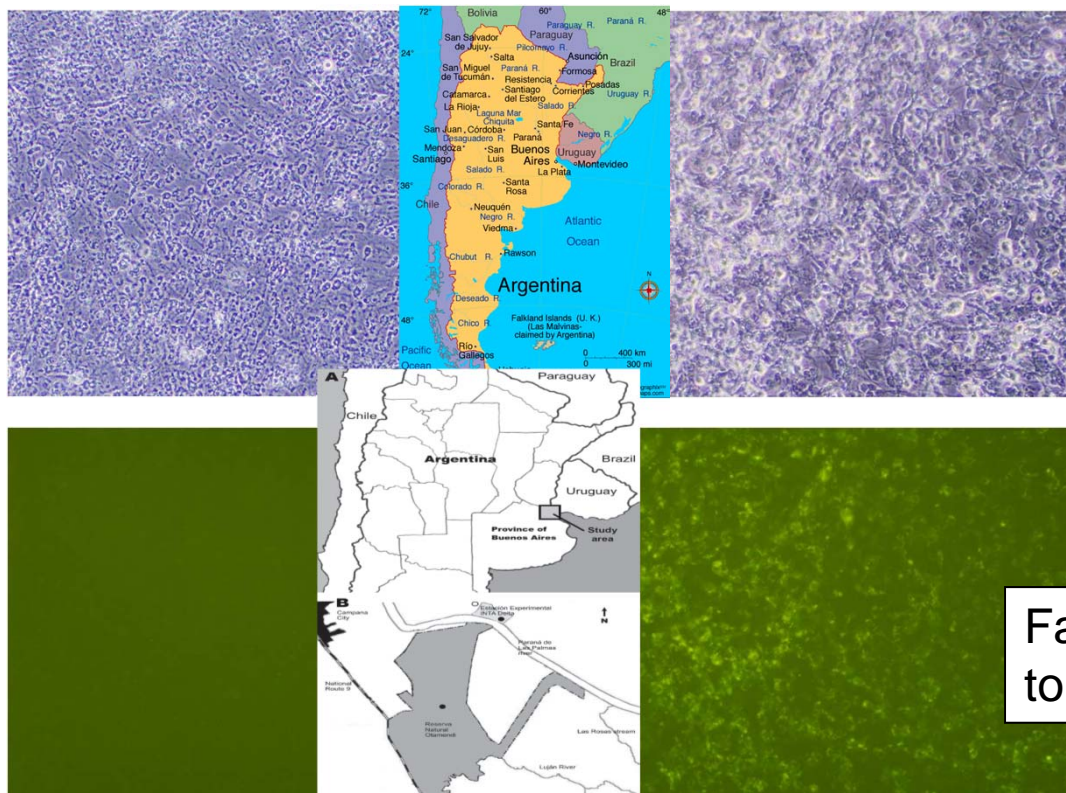

Suspected RLO attached by OX-19 aptamer to
nanocrystal of iron oxide on micro mag bead

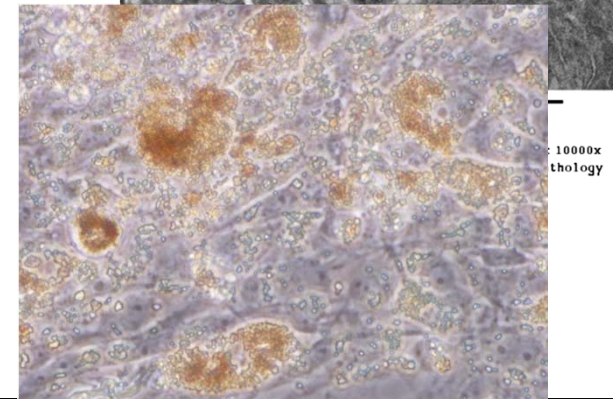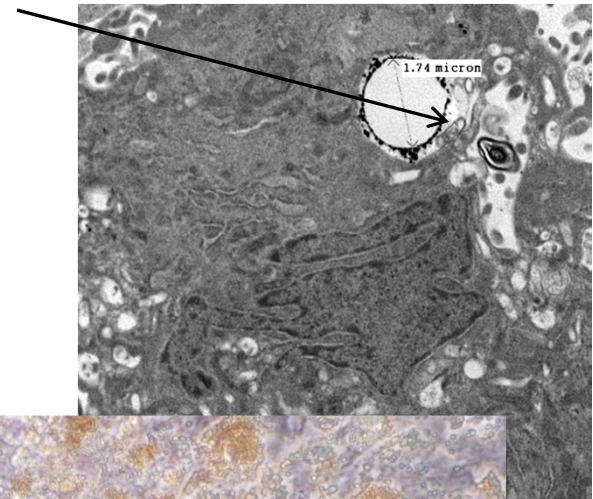

Facilitated Uptake of RLO Bound
to Beads (to bursting of cells)

I. Rosas and D. Vela

J. Kiel, R. Alarcon, J. Parker, J. Vivekananda, Y. Gonzalez, L. Stribling, C. Andrews, ANYAS **1081**: 434-442 (2006)
Kiel, J.L., Gonzalez, Y., Parker, J.E., Andrews, C., Martinez, D., Vacheiry, N., LeFrancois, T., ANYAS **1149**, 318-321 (2008).

Approved for public release; distribution unlimited

Summary

- Aptamers need to be selected under the conditions in which they are going to be used
 - **SELEX** aptamers sometime work as double-stranded contact reporting aptamers, but many times do not in spite of very low Kds
 - **ASExpP** fulfills the above criteria
- **SELEX**, by its very nature and mass action, selects for aptamers against the most abundant ligand not necessarily the most specific
 - **ASExpP**, because of its low cycle number and initial stringent conditions, selects for the highest affinity aptamer to the rarest target
- Several photochemical and electronic options exist for sensing platforms for aptamers
- Rapidity of aptamer selection in general allows for fast response to new emerging agents
- Double-stranded DNA capture elements allow for detection, identification and non-destructive safe collection for further orthogonal analysis in the lab
- Could emerging “Out of Africa” type D retrovirus be co-infecting with rickettsia transmitted by a tick vector and rendering the host more susceptible to the rickettsia by immunosuppression?
 - **A mechanism for emerging infectious disease ?**
 - **This new approach to environmental and animal sample collection may answer the question**

Approved for public release; distribution unlimited